2026 Product Presentation

LHALA JET

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

FRACTIONAL PLASMA **Low-temperature Atmospheric Pressure**

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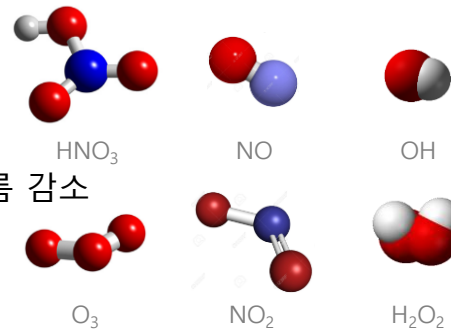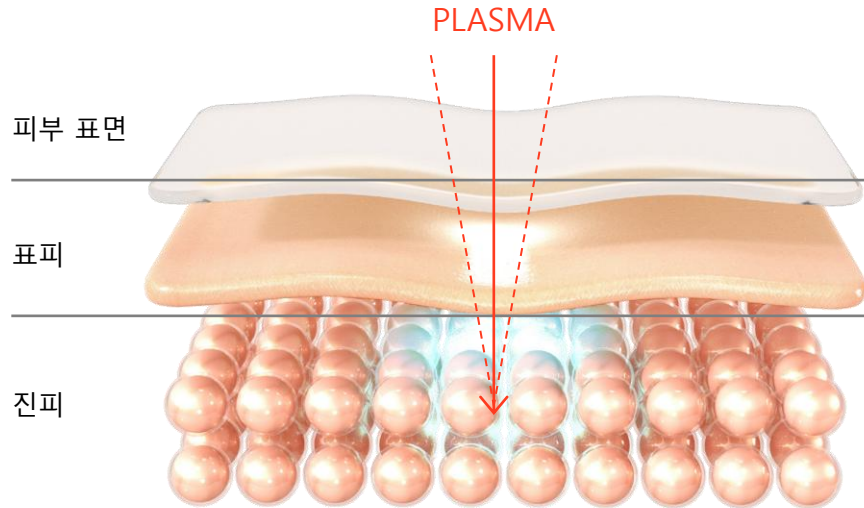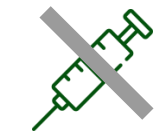
Drug Metabolism and Pharmacokinetics 54 (2024) 100536

Drug Metabolism and Pharmacokinetics

journal homepage: www.journals.elsevier.com/drug-metabolism-and-pharmacokinetics

Regular Article



Improvement of transdermal absorption rate by nonthermal biocompatible atmospheric pressure plasma

Byoung-Choul Kim^{a,1}, Juie Nahushkumar Rana^{a,b,1}, Eun Ha Choi^{a,b,c,**}, Ihn Han^{a,b,*}

ABSTRACT

Nonthermal biocompatible plasma (NBP) is a promising option for improving medication absorption into the human skin. Currently, most plasma devices for cosmetics employ a floating-electrode plasma source for treating the skin. Human skin serves as the ground electrode in the floating-electrode plasma discharge, and discharge occurs between the skin and electrodes of the device. In this *in vitro* study, we aimed to evaluate the effect of NBP on the skin permeation of niacinamide. We have quantified the transdermal absorption rates of niacinamide in both untreated skin and skin treated with NBP for a duration of 10 s. The absorption of niacinamide for both without and with NBP treatment was observed until 12 h incubation time. Without plasma treatment, the human skin exhibited stable and low transdermal absorption of niacinamide up to 12 h. However, the NBP treatment significantly increased the transdermal absorption of niacinamide from 0.5 h to 6 h and continuously increased skin penetration over a duration of more than 12 h incubation period. The obtained results suggest that **NBP-treated human skin showed a 60-fold higher penetration rate than non-treated skin**. The increased penetration rate of niacinamide can be mainly attributed to plasmaporation subsequent to NBP treatment. The findings of this study demonstrate that NBP treatment results in remarkable skin permeability, making it a promising candidate for both cosmetic and pharmaceutical delivery applications.

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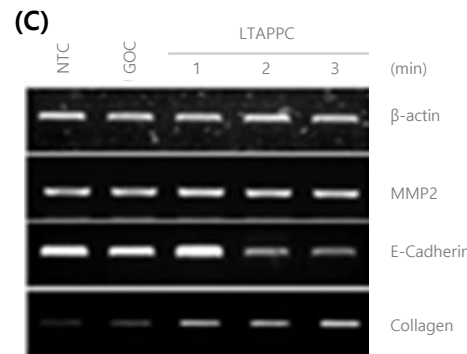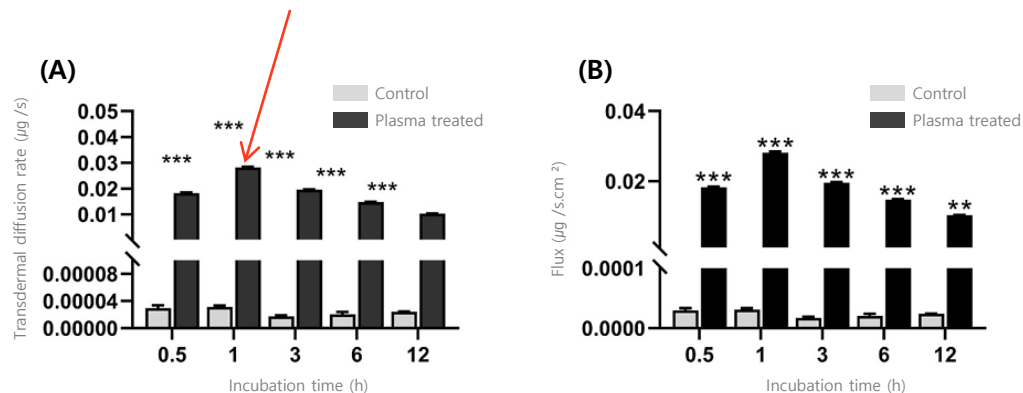Research Paper  Check for updates

Clinical Study on the Application of Non-Thermal Biocompatible Plasma for Acne-Prone Skin Treatment

Received 27 June, 2024; revised 2 July, 2024; accepted 2 July, 2024

Byoung Choul Kim^a, Youngsun Kim^b, Eun Ha Choi^{a,c,*}, and Ihn Han^{a,c,*}

^aDepartment of Plasma Bio Display, Kwangwoon University, Seoul 01897, Republic of Korea
^bDepartment of Obstetrics and Gynecology, Kyung Hee University Medical Center, Seoul 02447, Republic of Korea
^cPlasma Bioscience Research Center, Kwangwoon University, Seoul 01897, Republic of Korea

*Corresponding author E-mail: ehchoi@kw.ac.kr, hanihn@kw.ac.kr

ABSTRACT

Non-thermal biocompatible plasma (NBP), which generates various reactive species, is being actively researched in the field of aesthetic medicine for applications such as microorganism inactivation, skin regeneration, inflammation treatment, and drug absorption. This study evaluates the effectiveness of NBP and cosmetics in treating acne-prone skin by utilizing the sterilizing and regenerative effects of plasma. The present study involves 21 adult men and women over the age of 20 who require acne treatment. Participants use NBP and cosmetics after washing their faces in the morning and evening for six weeks. After six weeks, the global acne grading system scores show a 26.91 % reduction in acne severity and a 25.11 % reduction in sebum levels. The results confirm that NBP is effective in improving acne and reducing sebum. No adverse skin reactions are observed during the study period. Notably, NBP applies floating electrode-dielectric barrier discharge technology, which promotes the removal of acne-causing bacteria and skin regeneration. This study suggests that NBP is an effective method for acne treatment.



Keywords: Non-thermal biocompatible plasma, Floating electrode-dielectric barrier discharge, Acne treatment, Skin rejuvenation

Table I. Changes in acne grade score and improvement rate (%) after two, three, and six weeks of use (N=21). The results are presented as the mean ± standard deviation (SD) of three individual experiments. **p* < 0.05, ***p* < 0.01, and ****p* < 0.001.

	GAGS score			
	Control	2 weeks	3 weeks	6 weeks
Average	10.62	9.24	8.29	7.76
SD	3.31	3.16	2.53	2.72
Improvement rate		13.00	21.97	26.91
<i>p</i> -value		0.030*	0.001**	0.000***

Table II. Changes in sebum level and sebum improvement rate (%) after two, three, and six weeks of use (N=21). The results are presented as the mean ± SD of three individual experiments. ***p* < 0.01 and ****p* < 0.001.

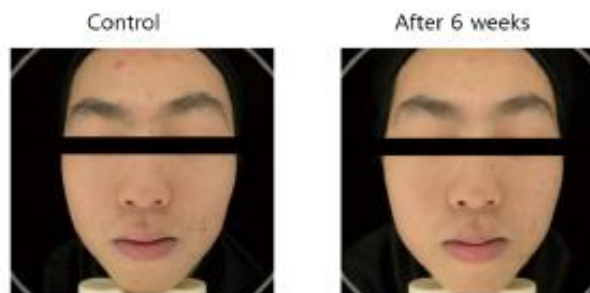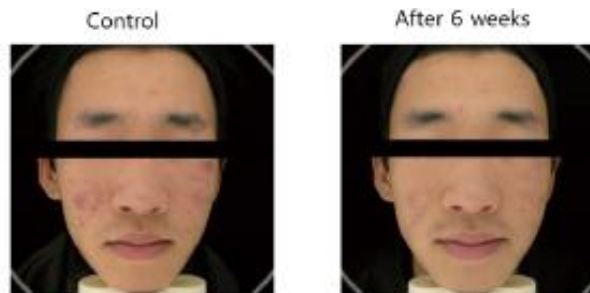	Control	Two weeks	Three weeks	Six weeks
Average	12.90	11.95	11.33	9.66
SD	6.25	6.25	6.19	6.02
Improvement rate		7.35	12.15	25.11
<i>p</i> -value		0.001**	0.000***	0.000***

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